

## 07 Series Fridge Controller Guide

This is the power button on the controller:



This button is to turn the light on. The interior light will stay on for two minutes or until you operate the button again to turn it off.



This is the key lock button and when it is lit up all other buttons other than the power button are inoperable.



This display shows the temperature. This reading is showing 5°.



This button is to check the low temperature. Press the button to display the lowest recorded temperature since the controller was last reset. This low temperature is reading -2°




This button is to check the high temperature. Press the button to display the highest recorded temperature since the controller was last reset. This high temperature is reading 5°



To reset the recorded high and low temperatures in memory, **press and hold** either of the temperature buttons and then press and hold the other temperature button simultaneously for 3 seconds until you hear a 2<sup>nd</sup> beep.



To reset the recorded high and low temperatures in memory, **press and hold** either of the temperature buttons and then press and hold the other temperature button simultaneously for 3 seconds until you hear a 2<sup>nd</sup> beep.

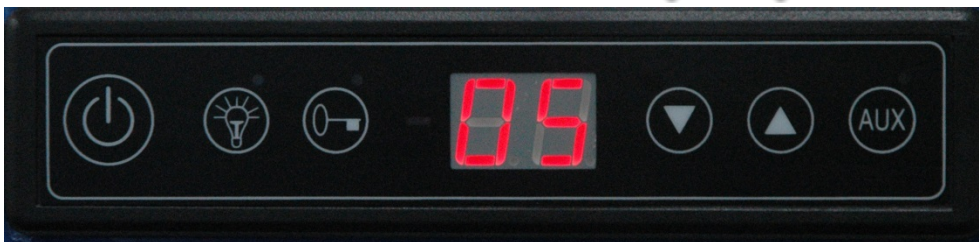

When the display reads “dr” it means that the door is open. The controller will start beeping if the door is open for one minute or longer and will stop when the door is closed or the door switch is depressed.



When the display reads “HI” and the light comes on this means the high alarm is on. The audible alarm will sound for approximately 3 minutes when an alarm condition is breached, the display will continue to show the alarm condition and the aux lamp will stay lit until a button on the controller is depressed.



If the high alarm is going off it may be from the door being open or a power cut if one has occurred. If it starts beeping again after you have turned the alarm off, please get in contact with us.

**Rollex Medical Head Office**

Phone: 0800 249 3722    Email: [nzservice@rollexmedical.co.nz](mailto:nzservice@rollexmedical.co.nz)

If you would like any further clarification on any of these matters please get in touch with us.

**Rollex Medical NZ Limited**

Postal Address: PO Box 204159, Highbrook, Auckland 2161, New Zealand

Physical Address: 38E Highbrook Drive, East Tamaki, Auckland 2013, New Zealand

Freephone: 0800 249 3722

Phone: +64 9 274 0707

Fax: +64 9 274 0709

Email (Sales): [nzsales@rollexmedical.co.nz](mailto:nzsales@rollexmedical.co.nz)

Email (Service): [nzservice@rollexmedical.co.nz](mailto:nzservice@rollexmedical.co.nz)

